

OL 730-9A and OL 730-9Q Reflex Telescopes

The OL 730-9A and OL 730-9Q Reflex Telescopes enable the Optronic Laboratories series of light measuring instruments to operate either as telephotometers, teleradiometers, or telespectroradiometers. The reflex telescope is used to measure the luminance, radiance, or spectral radiance of uniform, diffusely emitting light sources.

OL 730-9A

The OL 730-9A consists of an achromatic, 66 mm lens in a rack and pinion focusing mount, a baffle tube, and a reflex viewer compartment with zero parallax viewing. The reflex viewer has a three element, orthoscopic eyepiece that yields an 8X magnification. The wide-field eyepiece has a precision reticle with ten concentric circles ranging from 1 to 10 mm. Three apertures producing FOVs (field of views) of 1.5°, 1°, and 0.5° are standard. The useful wavelength range of the OL 730-9A is 360 to 1100 nm.

OL 730-9Q

The OL 730-9Q is identical to the OL 730-9A except for the substitution of a quartz lens for the standard achromatic lens. The quartz lens extends the useful wavelength transmission range from beyond 250 nm in the UV to 2.5 μm in the IR. However, this is accomplished with a slight deterioration of the focused image due to chromatic aberration.


A mounting flange for attaching the reflex telescope directly to the monochromator and a removable adapter ring for coupling the detector (with the appropriate filter) directly to the telescope are supplied. Fittings for use with the OL 730-7 Fiber Optic Probe are also available.

Approximate photometric sensitivities when the OL730-9A/Q is used with the OL 730/730A/730C Radiometer/Photometer are:



FOV Aperture	Noise Equivalent Luminance (fL)	
	W/ silicon detector	W/ PMT detector
1.5°	3×10^{-4}	5×10^{-7}
1.0°	1×10^{-3}	2×10^{-6}
0.5°	4×10^{-3}	8×10^{-6}

Approximate spectroradiometric sensitivities when the OL 730-9A/Q, with the 1° FOV aperture, is used with the OL 740A, OL 746, or OL 750 Spectroradiometer are:

Wavelength	Noise Equivalent w/ Silicon Detector	(Watt/sr cm ² nm) w/ PMT Detector
360 nm	1×10^{-10}	8×10^{-14}
500 nm	5×10^{-11}	2×10^{-14}
600 nm	6×10^{-11}	5×10^{-14}
700 nm	7×10^{-11}	1×10^{-13}
800 nm	8×10^{-11}	5×10^{-13}
1000 nm	9×10^{-11}	N/A

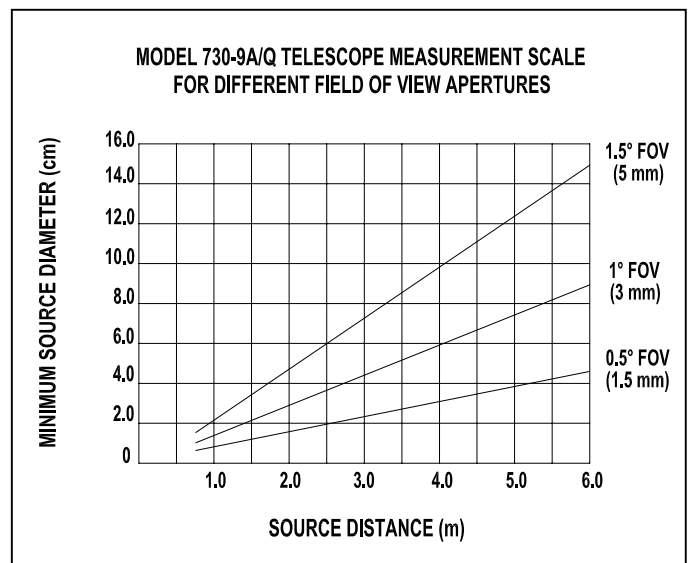

OL 730-9A and OL 730-9Q
 Reflex Telescopes

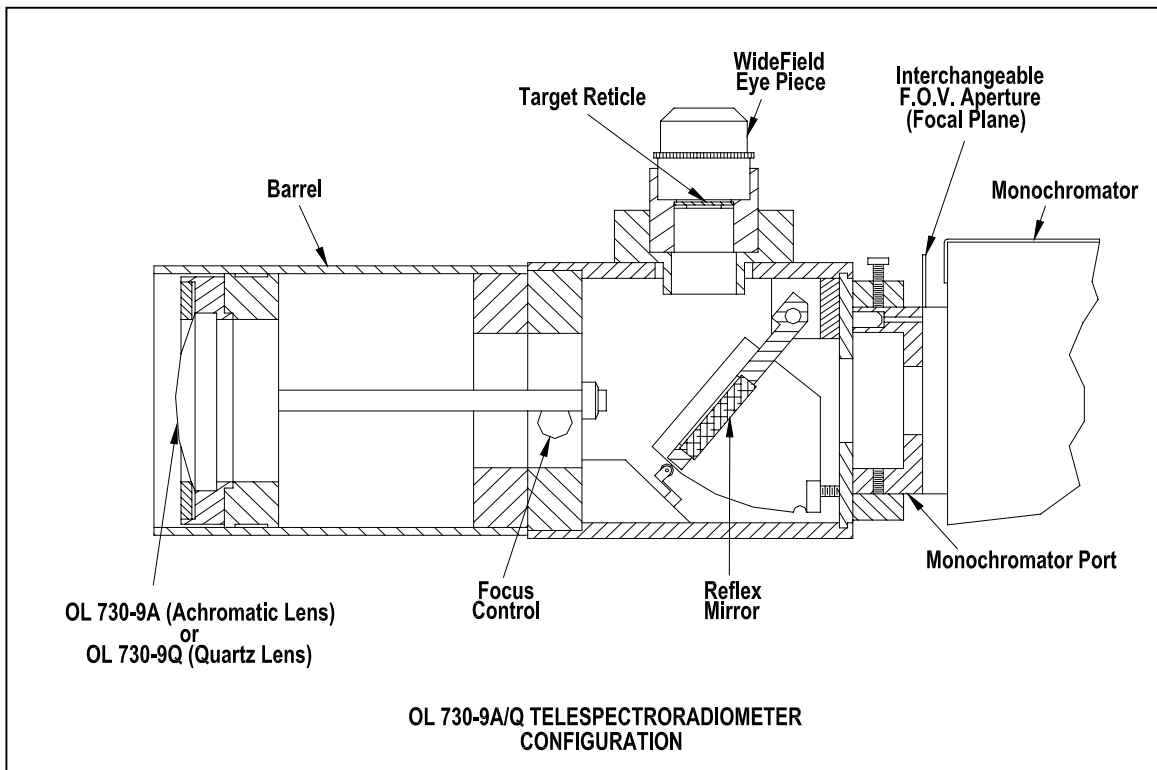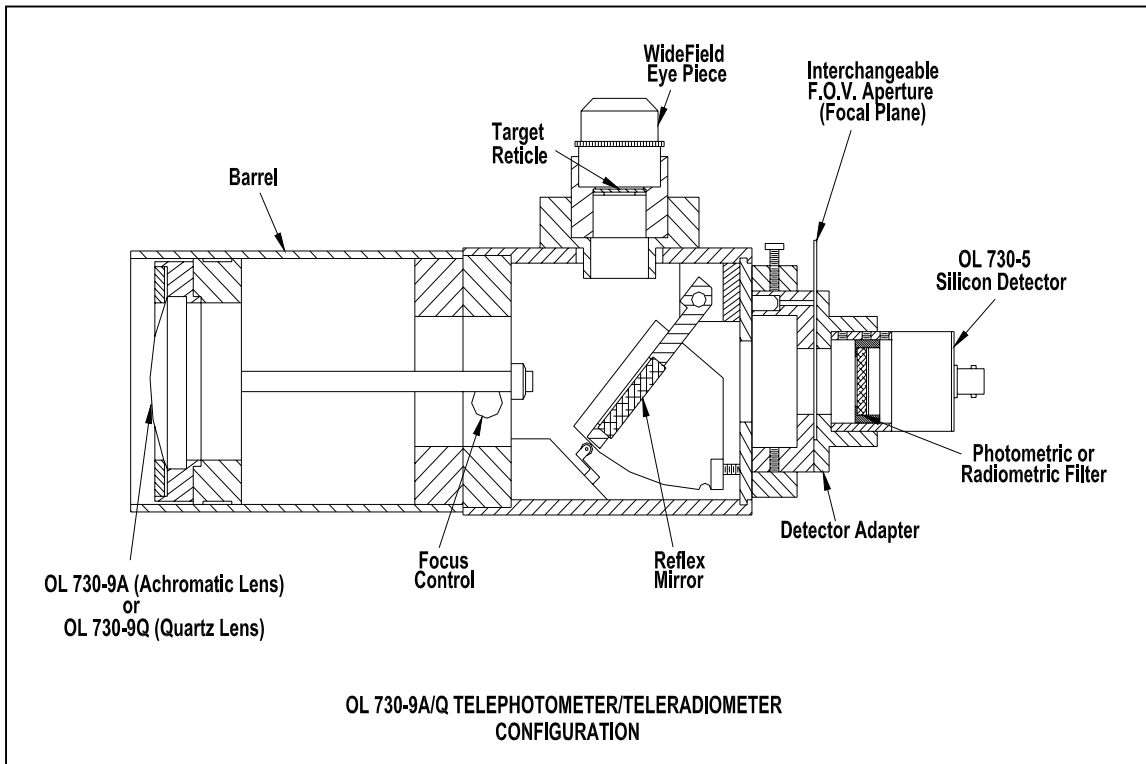

OPTRONIC LABORATORIES

A Gooch & Housego Company

Optronic Laboratories, Inc. 4632 36th Street, Orlando, FL 32811
 Tel: 1 407 422 3171 Fax: 1 407 648 5412 Email: info@olinet.com

Bulletin 115 / Rev. 3-02